

Effects of bonding materials on the mechanical properties of enamel around orthodontic brackets

Naohisa Kohda^a; Masahiro Iijima^b; William A. Brantley^c; Takeshi Muguruma^a; Toshihiro Yuasa^a; Susumu Nakagaki^d; Itaru Mizoguchi^e

ABSTRACT

Objective: To determine if the enamel around orthodontic brackets is significantly altered after demineralization followed by application of adhesives with and without fluoride-releasing ability.

Materials and Methods: One hundred eight noncarious human premolars were divided into six groups of 18 each and exposed to a demineralization solution. Stainless steel brackets were bonded using two conventional composite resin etch-and-rinse systems, three self-etching primer (SEP) composite resin systems, and one resin-modified glass ionomer cement (RMGIC) system. One conventional and one SEP composite resin adhesive did not have fluoride-releasing ability, which was claimed for the other four adhesives. The elastic modulus and hardness of the enamel were determined with a nanoindenter at 10 equidistant depths ranging from 1–46 μm and at four regions: control (not exposed) enamel surface, under the adhesive, and at 50 μm and 100 μm from the bracket edges. Using the Kruskal-Wallis and Mann-Whitney *U*-tests ($P < .0125$ for statistical significance), these properties were compared at different regions.

Results: The same behavior was observed for values of elastic modulus and hardness. Significant differences were found within approximately 21 μm of the enamel surface for etching with 35% phosphoric acid or priming with SEP, but only minimal changes occurred for the SEP adhesive. Increases in near-surface elastic modulus and hardness of enamel were found with the SEP adhesive and RMGIC with fluoride-releasing ability.

Conclusions: Clinical use of the fluoride-releasing adhesives investigated may prevent demineralization of enamel around brackets during orthodontic treatment. (*Angle Orthod.* 2012;82:187–195.)

KEY WORDS: Enamel; Demineralization; Nanoindentation; Hardness; Elastic modulus; Fluoride-releasing adhesive

INTRODUCTION

Increased prevalence of white spot formation in orthodontic patients has been reported¹ due to the irregular surface of orthodontic fixed appliances, which creates stagnation areas for plaque, render tooth cleaning more difficult, and limits naturally occurring self-cleaning mechanisms.¹ Fluoride plays an important role in the prevention of demineralization during orthodontic treatment, and mouth rinses with fluoride solutions are recommended in addition to daily tooth-brushing. However, the use of mouth rinses for preventive enamel demineralization requires cooperation by the patient, which is sometimes difficult, especially for young patients. On the other hand, the original method of using phosphoric acid for etching in bracket bonding is

^a Instructor, Division of Orthodontics and Dentofacial Orthopedics, Department of Oral Growth and Development, Health Sciences University of Hokkaido, Hokkaido, Japan.

^b Associate Professor, Division of Orthodontics and Dentofacial Orthopedics, Department of Oral Growth and Development, Health Sciences University of Hokkaido, Hokkaido, Japan.

^c Professor, Division of Restorative and Prosthetic Dentistry, College of Dentistry, The Ohio State University, Columbus, Ohio.

^d PhD student, Division of Orthodontics and Dentofacial Orthopedics, Department of Oral Growth and Development, Health Sciences University of Hokkaido, Hokkaido, Japan.

^e Professor and Chairman, Division of Orthodontics and Dentofacial Orthopedics, Department of Oral Growth and Development, Health Sciences University of Hokkaido, Hokkaido, Japan.

Corresponding author: Dr Masahiro Iijima, Division of Orthodontics and Dentofacial Orthopedics, Department of Oral Growth and Development, Health Sciences University of Hokkaido, 1757, Kanazawa, Ishikari-tobetsu, Hokkaido 061-0293 Japan (e-mail: iijima@hoku-iryō-u.ac.jp)

Accepted: June 2011. Submitted: January 2011.

Published Online: August 9, 2011

© 2012 by The EH Angle Education and Research Foundation, Inc.

associated with a loss of the enamel surface (approximately 10–20 μm),^{2,3} although etched enamel around orthodontic brackets might be remineralized by an application of the fluoride products.⁴

Over the past decade, progress has been made in bracket bonding to enamel with resin-modified glass ionomer cements (RMGICs) and self-etching primers (SEPs),² and their lower etching ability might minimize the potential for iatrogenic damage to enamel. Fluoride-releasing adhesives with less etching ability, such as RMGIC and SEP composite resin adhesive systems containing fillers with fluoride, might be good selections for preventive enamel demineralization. The preventive effect of an adhesive adjacent to brackets has been investigated *in vitro*^{1,5,6} and *in vivo*^{7–9} by quantifying the demineralization depths and the mineral losses with various evaluation methods. Most studies have been shown that RMGICs are more effective for prevention of enamel demineralization in comparison with fluoride-releasing and non-fluoride-releasing composite resin adhesive systems.^{1,5,8} To investigate the demineralization of enamel quantitatively, use of cross-sectional microhardness measurements with a Knoop indenter has been a popular method^{5–8} because there is high correlation between enamel microhardness and the percentage of mineral in carious lesions.¹⁰

Recent advances in the nanoindentation test have allowed the measurement of mechanical properties for extremely small volumes of materials where the contact radius is less than 100 nm.^{11–14} In nanoindentation testing, the load and displacement are recorded with a high-resolution displacement gage during the indentation process, and the area of the indentation is obtained from the geometric form of the indenter without having to actually observe the indentation; thereby the hardness of the specimen can be calculated. In addition, the elastic modulus for this very small volume of material can be obtained mathematically from the load-displacement curve.

Little is presently known about mechanical properties of the subsurface region of enamel around orthodontic brackets. The purpose of this *in vitro* study was to investigate the hardness and elastic modulus of enamel around orthodontic brackets after immersion in a demineralization solution and compare these mechanical properties for different adhesive systems with and without fluoride-releasing ability. It was hypothesized that different types of adhesive systems have significantly different effects on the enamel demineralization behavior and mechanical properties.

MATERIALS AND METHODS

Materials

One hundred eight noncarious human premolars, obtained by extraction from patients who were to

undergo orthodontic treatment, were allocated into six groups of 18 each. This *in vitro* study received approval by an ethics committee at the Health Sciences University of Hokkaido. The buccal surfaces of all teeth were cleaned using nonfluoridated pumice. The teeth were subsequently polished using a rubber cup, and thoroughly washed and dried using a moisture-free air source. Maxillary or mandibular stainless steel brackets (Victory Series, 3M Unitek, Monrovia, Calif) were bonded to the enamel samples with each of the six bracket bonding adhesive systems listed in Table 1. Groups 1 and 2 employed conventional etch-and-rinse composite resin adhesive systems. Groups 3, 4, and 5 employed SEP composite resin adhesive systems. Group 6 employed an RMGIC adhesive system. The manufacturers claim that the adhesives employed for groups 2, 4, 5, and 6 have fluoride-releasing ability.

Preparation of Teeth

The buccal surfaces of the premolars were coated with a 1.5-mm width of acid-resistant nail varnish from the cusp to the cervix, and these regions under the nail varnish were used as the control (Figure 1a). In groups 1 and 2, the buccal surfaces were etched with conventional 35% phosphoric acid gel for 15 seconds (Figure 1b), washed for 20 seconds, and dried with an oil-free air stream. Transbond XT primer was then applied to the etched surface. In groups 3, 4, and 5, Transbond SEP or Beauty Ortho Bond SEP was applied by rubbing on the buccal enamel surfaces for approximately 3 seconds; afterward an air jet was lightly applied. In group 6, 20% polyacrylic acid Ortho Conditioner was applied for 10 seconds on the buccal enamel surface, which was washed for 20 seconds, and dried with an oil-free air stream. Then, the stainless steel brackets (Victory Series, 3M Unitek) were bonded with the composite resin adhesive pastes or RMGIC paste and light-cured according to the instructions from each manufacturer (Figure 1c). Acid-resistant nail varnish was next applied to each tooth, leaving a 1-mm rim of exposed sound enamel surrounding the bracket (Figure 1d). The specimens were then immersed in a demineralization solution (4% Methocel MC gel [Fluka, Buchs, Switzerland] + 0.1% lactic acid) in a polystyrene vial (Figure 1e). The solution was adjusted to a pH of 4.6. All specimens were immersed in 10 mL of this solution for 96 hours at 37°C, with the solution being changed every 12 hours, and then gently washed with deionized water.

All of the teeth were cut with a slow-speed water-cooled diamond saw (Isomet, Buehler, Lake Bluff, Ill) so that they were divided into occlusal and cervical halves (Figure 1f); one of the sectioned specimens

Table 1. Materials Employed in Present Study

	Material ^a	Components	Composition	Fluoride	Manufacturer
Group 1	Transbond XT etch-and-rinse adhesive	Etching gel Primer Paste	35% Phosphoric acid TEGDMA, Bis-GMA TEGDMA, Bis-GMA, silane-treated quartz, amorphous silica, camphorquinone		3M Unitek
Group 2	Transbond Plus Color Change etch-and-rinse adhesive	Etching gel Primer Paste	35% Phosphoric acid TEGDMA, Bis-GMA TEGDMA, Bis-GMA, silane-treated quartz, amorphous silica, camphorquinone	+	3M Unitek
Group 3	Transbond XT SEP adhesive	SEP Paste	Methacrylated phosphoric acid esters, amino benzoate, camphorquinone TEGDMA, Bis-GMA, silane-treated quartz, amorphous silica, camphorquinone		3M Unitek
Group 4	Transbond Plus Color Change SEP adhesive	SEP Paste	Methacrylated phosphoric acid esters, amino benzoate, camphorquinone TEGDMA, Bis-GMA, silane-treated quartz, amorphous silica, camphorquinone	+	3M Unitek
Group 5	Beauty Ortho Bond SEP adhesive	Primer A Primer B Paste	Water, acetone, other Phosphoric acid monomer, ethanol TEGDMA, Bis-GMA, S-PRG filler, camphorquinone	+	Shofu
Group 6	Fuji Ortho LC Automix RMGIC adhesive	Conditioner Paste A Paste B	20% Polyacrylic acid Fluoro-alumino-silicate glass, dimethacrylate polyacrylic acid, silicon dioxide, distilled water	+	GC America

^a SEP indicates self-etching primer; RMGIC, resin-modified glass ionomer cement.

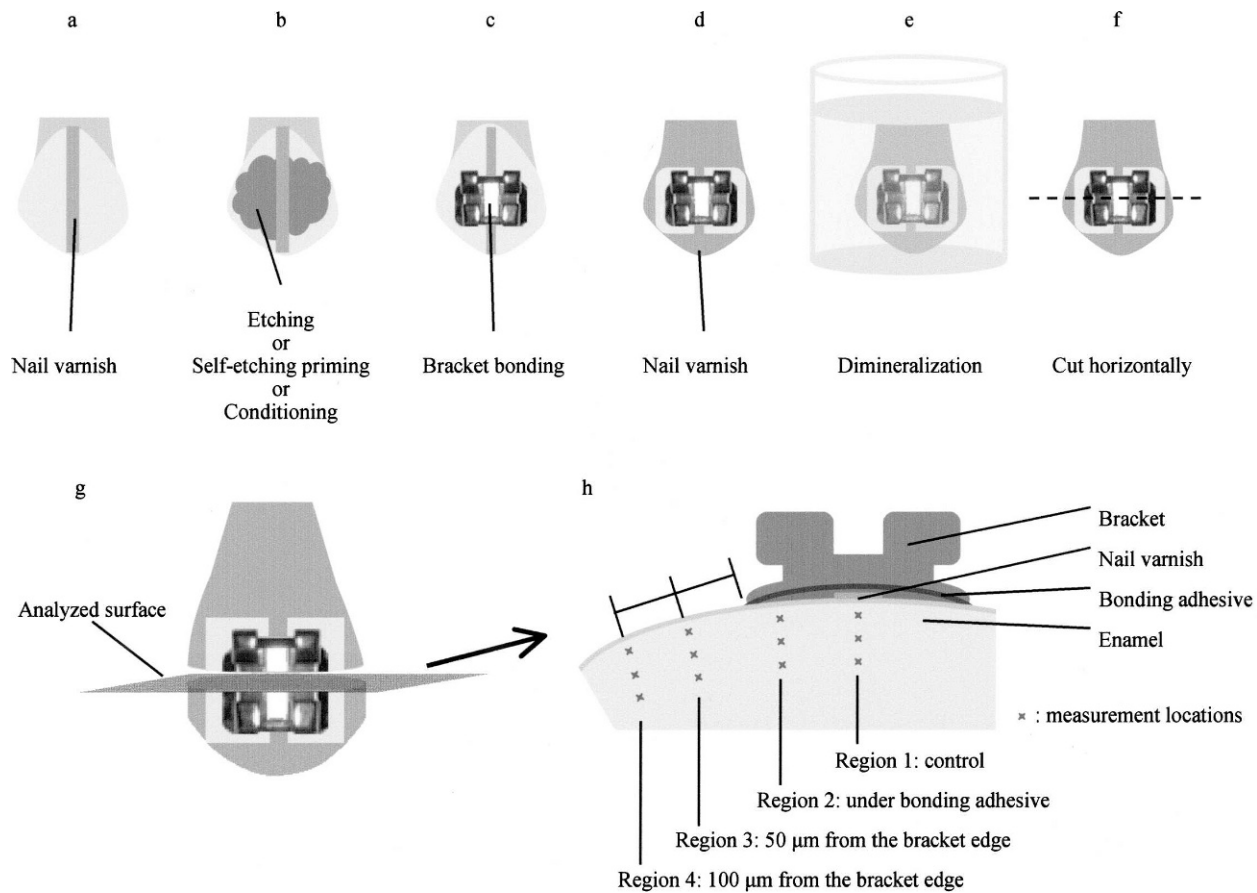


Figure 1. Schematic illustration of specimen preparation sequence and regions of tooth enamel investigated with the nanoindentation test.

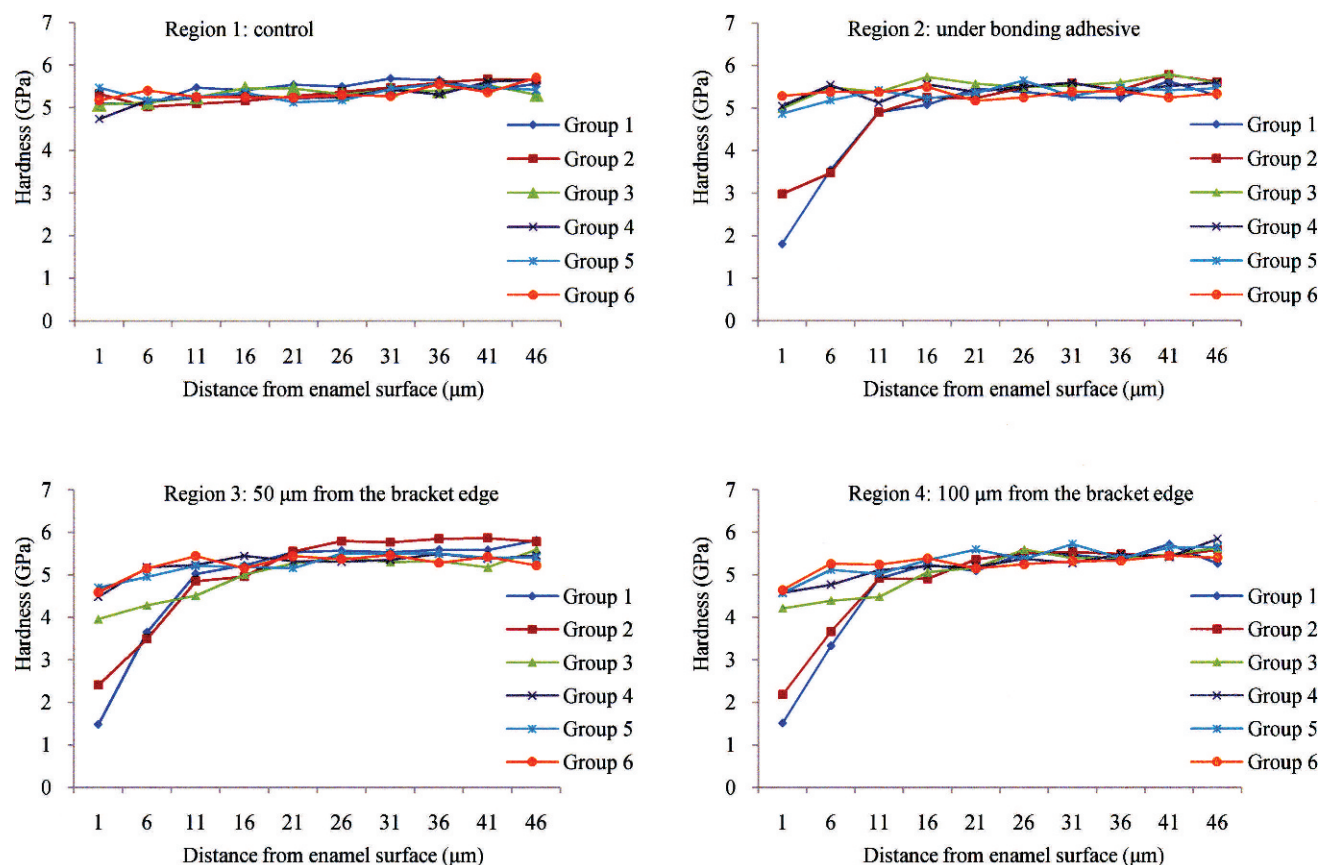


Figure 2. Mean values of hardness of enamel at different distances from the surface after immersion in demineralization solution. Group 1, Transbond XT etch-and-rinse adhesive system; group 2, Transbond Plus Color Change etch-and-rinse adhesive system; group 3, Transbond XT SEP adhesive system; group 4, Transbond Plus Color Change SEP adhesive system; group 5, Beauty Ortho Bond SEP adhesive system; and group 6, Fuji Ortho LC Automix RMGIC adhesive system.

(transverse planes) was then encapsulated in epoxy resin (Epofix, Struers, Copenhagen, Denmark) for the nanoindentation test. All samples were ground (600-grit sandpaper) and polished using diamond suspensions (3-, 1-, and 0.25-μm particle sizes) to obtain a suitable surface for nanoindentation.

Nanoindentation Test

All nanoindentation testing (ENT-1100a, Elionix, Tokyo, Japan) was carried out at 28°C with a peak load of 10 mN using a Berkovich indenter. Each test consisted of three parts: 10 seconds for loading to the peak value, 1 second of holding at the peak load, and 10 seconds for unloading. The indentations were placed at 1–46 μm depth (10 locations spaced 5 μm apart) from the external surface for four regions (region 1: control; region 2: under the bonding adhesive; region 3: 50 μm from the bracket edge; region 4: 100 μm from the bracket edge) (Figure 1h). The measurement points were selected using the optical microscope and charge-coupled device camera coupled to the nanoindentation test apparatus. Linear extrapolation methods (ISO Standard 14577-1)¹⁵ were used for the unloading curve

between 95% and 70% of the maximum test force to calculate the elastic modulus. The hardness and elastic modulus were calculated by software available with the nanoindentation apparatus.

Statistical Analysis

The experimental results were analyzed using the PASW Statistics software (version 18.0J for Windows, IBM, Armonk, New York). The data for the mean values of hardness and elastic modulus were not normally distributed (Levene test), so the values were compared using the Kruskal-Wallis test and the Mann-Whitney *U*-test ($P < .0125$). The purpose of these comparisons was to assess the effect of bracket bonding procedure or immersion in demineralization solution (regions 1–4) for each adhesive.

RESULTS

Figures 2 and 3 show mean values for hardness and elastic modulus of the cross-sectioned polished specimens (transverse planes) obtained by the nanoindentation test; the results of statistical comparisons

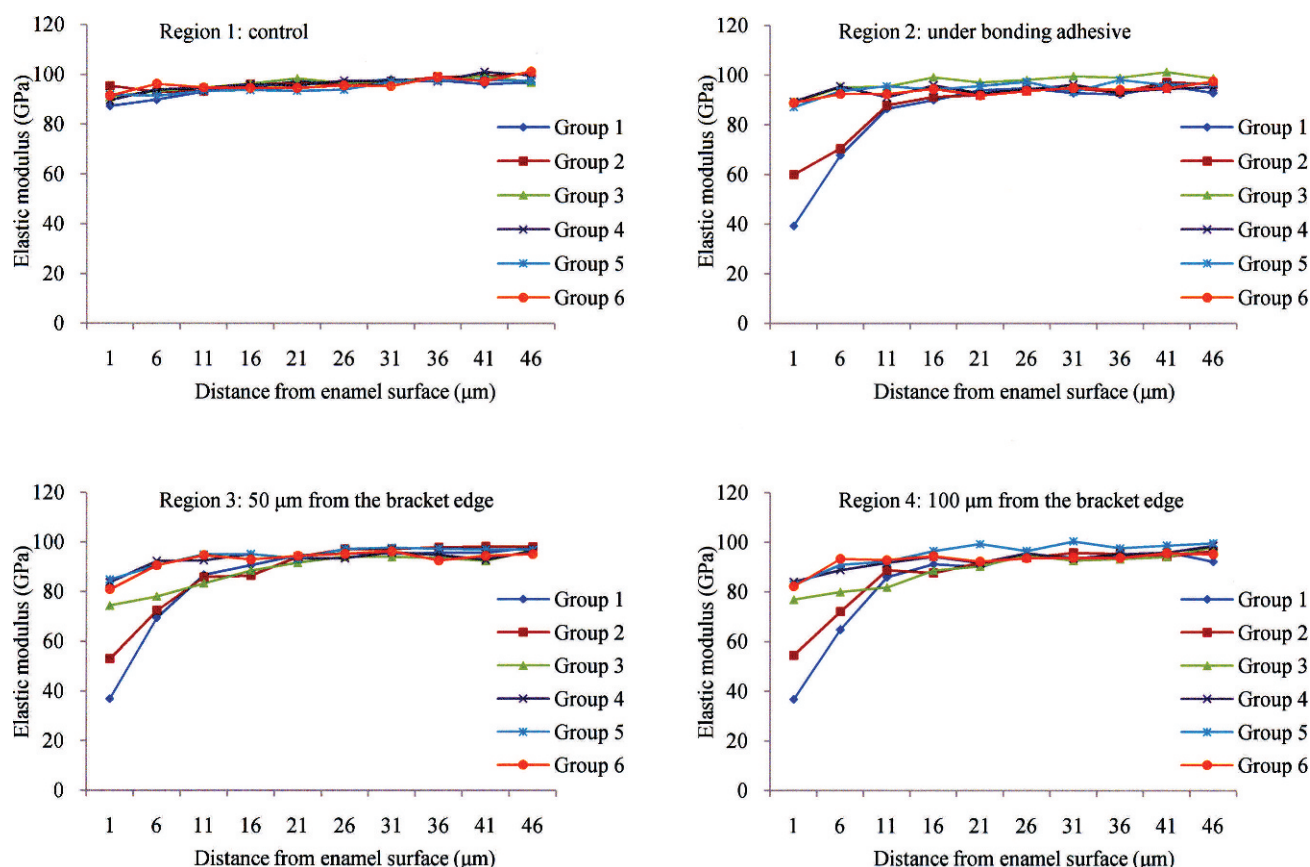


Figure 3. Mean values of elastic modulus of enamel at different distances from the surface after immersion in demineralization solution. Group 1, Transbond XT etch-and-rinse adhesive system; group 2, Transbond Plus Color Change etch-and-rinse adhesive system; group 3, Transbond XT SEP adhesive system; group 4, Transbond Plus Color Change SEP adhesive system; group 5, Beauty Ortho Bond SEP adhesive system; and group 6, Fuji Ortho LC Automix RMGIC adhesive system.

between different regions are summarized in Tables 2 and 3. Locations at 1 μm and 6 μm from the enamel surface for regions 2, 3, and 4 of both conventional etch-and-rinse adhesive systems (groups 1 and 2) had significantly lower hardness and elastic modulus than those for region 1, except for the hardness of region 4 in group 2 at 6 μm depth. For group 3, there was no significant difference in the hardness and elastic modulus between regions 1 and 2. The hardness up to 16 μm in depth and elastic modulus up to 21 μm in depth from the enamel surface for regions 3 and 4 had significantly lower values than those for region 1 and 2 except for elastic modulus at 16 μm depth from the enamel surface for region 4. There was no significant difference in the mean values for hardness and elastic modulus of enamel at each depth for group 4. For group 5, the locations at 1 μm depth from the enamel surface for regions 2, 3, and 4 had significantly lower hardness than region 1, although there were no significant differences in elastic modulus. For group 6, the location at 1 μm below the enamel surface for region 3 had significantly lower hardness than region 2, but region 2 was not significantly different from region 1, and region

3 was not significantly different from region 4. The elastic modulus for group 6 at 1 μm depth was significantly different for regions 1 and 3, but not for regions 1 and 2 nor for regions 3 and 4.

DISCUSSION

The present in vitro study evaluated the effect of bonding materials on the mechanical properties (hardness and elastic modulus) of enamel around orthodontic brackets. These mechanical properties may change depending upon the bracket bonding procedure and demineralization of enamel, and in this study changes were assessed by cross-sectional nanoindentation analysis. The mechanical properties of enamel are generally linked to its mineral content. Knowledge of the mechanical properties of enamel such as the hardness and elastic modulus may lead to an improved understanding of the bracket bonding and demineralization behavior, along with the effects of new bonding materials. In the nanoindentation test, the minimum distance between indentations should be at least five times the largest indentation diameter.¹⁵ This study used a 10 mN peak load, which made

Table 2. Mean Values and Standard Deviation for Cross-sectional Hardness (GPa) of Enamel at Distances From the Surface^a

Depth, μm	Region 1		Region 2		Region 3		Region 4	
	Mean	SD	Mean	SD	Mean	SD	Mean	SD
Group 1 : Transbond XT								
1	5.10 ^a	0.92	1.81 ^b	1.30	1.48 ^b	0.63	1.51 ^b	0.90
6	5.09 ^a	0.94	3.55 ^b	1.14	3.65 ^b	1.16	3.33 ^b	1.60
11	5.48 ^a	0.44	4.90 ^b	0.64	5.02 ^b	0.54	4.92 ^{a,b}	0.91
16	5.41	0.54	5.08	0.60	5.22	0.79	5.24	0.52
21	5.54	0.75	5.46	0.57	5.53	0.73	5.09	0.60
26	5.49	0.78	5.38	0.52	5.57	0.61	5.51	0.42
31	5.69	0.76	5.25	0.32	5.53	0.67	5.47	0.55
36	5.65	0.70	5.24	0.36	5.59	0.63	5.36	0.79
41	5.42	0.64	5.63	0.49	5.59	0.69	5.72	0.74
46	5.56	1.02	5.29	0.38	5.81	0.87	5.26	0.42
Group 2 : Transbond Plus Color Change								
1	5.33 ^a	0.81	2.98 ^b	1.56	2.41 ^b	1.15	2.18 ^b	1.32
6	5.02 ^a	0.73	3.48 ^b	1.39	3.49 ^b	1.40	3.66 ^{a,b}	1.60
11	5.10	0.66	4.90	0.64	4.85	0.72	4.91	0.75
16	5.16	0.73	5.25	0.39	4.96	0.62	4.90	0.83
21	5.27	0.87	5.23	0.53	5.55	0.47	5.36	0.77
26	5.37	0.58	5.49	0.42	5.79	0.53	5.51	0.59
31	5.47	0.53	5.58	0.58	5.77	0.35	5.53	0.71
36	5.59	0.53	5.42	0.48	5.85	0.34	5.49	0.91
41	5.67	0.87	5.78	0.68	5.86	0.33	5.44	1.15
46	5.65	0.63	5.61	0.75	5.78	0.31	5.59	0.80
Group 3 : Transbond XT + SEP								
1	5.08 ^a	0.34	5.00 ^a	0.56	3.96 ^b	0.95	4.21 ^b	0.95
6	5.12 ^a	0.40	5.49 ^a	0.46	4.28 ^b	0.84	4.39 ^b	0.91
11	5.25 ^a	0.38	5.38 ^a	0.55	4.51 ^b	0.90	4.48 ^b	0.98
16	5.46 ^{a,b}	0.38	5.73 ^a	0.53	5.01 ^b	0.97	5.06 ^b	0.80
21	5.45	0.56	5.57	0.44	5.27	0.69	5.17	0.57
26	5.32	0.51	5.49	0.67	5.41	0.66	5.59	0.56
31	5.41	0.55	5.53	0.63	5.30	0.58	5.39	0.62
36	5.40	0.52	5.60	0.45	5.32	0.82	5.40	0.59
41	5.53	0.60	5.79	0.78	5.18	0.61	5.47	0.68
46	5.30	0.69	5.58	0.64	5.59	0.65	5.62	0.80
Group 4 : Transbond Plus Color Change + SEP								
1	4.74	0.55	5.05	0.79	4.48	0.64	4.58	0.57
6	5.17	0.96	5.55	0.97	5.18	0.80	4.77	0.58
11	5.24	0.75	5.13	0.73	5.23	0.78	5.11	0.58
16	5.30	0.42	5.56	0.94	5.45	0.79	5.20	0.89
21	5.22	0.57	5.38	0.65	5.32	0.59	5.19	0.62
26	5.25	0.53	5.50	0.90	5.31	0.57	5.38	0.58
31	5.42	0.94	5.60	0.63	5.36	0.50	5.27	0.51
36	5.31	0.70	5.38	0.59	5.50	0.72	5.48	0.61
41	5.61	0.76	5.51	0.49	5.38	0.45	5.42	0.53
46	5.66	0.83	5.59	0.80	5.46	0.61	5.84	0.76
Group 5 : Beauty Ortho Bond								
1	5.47 ^a	0.54	4.87 ^b	0.53	4.70 ^b	0.85	4.57 ^b	0.82
6	5.17	0.71	5.18	0.66	4.95	0.68	5.12	0.80
11	5.26	0.74	5.42	0.78	5.21	0.64	5.02	0.56
16	5.36	0.70	5.23	0.65	5.19	0.72	5.34	0.89
21	5.13	0.69	5.35	0.66	5.16	0.70	5.59	0.65
26	5.18	0.68	5.65	0.71	5.50	0.88	5.37	0.52
31	5.44	0.59	5.28	0.48	5.51	0.52	5.72	0.87
36	5.57	0.61	5.48	0.52	5.50	0.78	5.39	0.50
41	5.47	0.49	5.42	0.41	5.40	0.45	5.64	0.54
46	5.42	0.36	5.46	0.57	5.41	0.38	5.65	0.98
Group 6 : Fuji Ortho LC Automix								
1	5.18 ^{a,b}	0.69	5.29 ^a	0.81	4.59 ^b	0.80	4.64 ^{a,b}	0.84
6	5.40	0.91	5.38	0.63	5.15	0.54	5.26	0.72

Table 2. Continued

Depth, μm	Region 1		Region 2		Region 3		Region 4	
	Mean	SD	Mean	SD	Mean	SD	Mean	SD
11	5.25	0.60	5.37	0.64	5.44	0.92	5.24	0.58
16	5.25	0.61	5.50	0.75	5.16	0.80	5.39	0.47
21	5.24	0.47	5.17	0.65	5.44	0.71	5.15	0.70
26	5.30	0.65	5.25	0.65	5.37	0.48	5.25	0.50
31	5.27	0.69	5.38	0.54	5.46	0.82	5.32	0.95
36	5.55	0.75	5.40	0.40	5.28	0.61	5.33	0.52
41	5.36	0.71	5.24	0.42	5.42	0.62	5.45	1.02
46	5.70	1.14	5.33	0.53	5.22	0.73	5.40	1.04

^a Region 1, control; region 2, under bonding adhesive; region 3, 50 μm from the bracket edge; region 4, 100 μm from the bracket edge. SEP indicates self-etching primer. Identical letters indicate that mean values were not significantly different ($P < .0125$).

indentations having about 2 μm width and length; the distance between the indentations was about 3 μm. Because a greater peak load shows more accurate values and stable data, and because the indentations obtained by our condition did not produce any cracks on the enamel, we used this 10 mN peak load condition. In addition, a distance of less than 5 μm apart between indentations was favorable to analyze the demineralization behavior in this study.

This in vitro study used noncarious human premolars obtained from different patients, and the observed individual variation in acid resistance of the enamel specimens would be expected. However, the hardness and elastic modulus for the control region (region 1) obtained from all adhesive systems (groups 1–6) showed no significant variation for all the locations except that the group 4 location at 1 μm had significantly lower hardness than the other five groups. Mean hardness and elastic modulus for region 1 ranged from 4.7–5.7 GPa and 87.3–99.6 GPa, respectively, and these values were similar to recently published values obtained with the nanoindentation test (3.4–8.3 GPa for hardness; 61–130 GPa for elastic modulus).^{3,12} The values of hardness and elastic modulus obtained in this study should be affected by the location (depth from enamel surface and analysis region for each adhesive group), along with test load, type of indenter, and direction of the enamel rods.

The locations at 1 μm and 6 μm from the enamel surface for region 2 for both conventional etch-and-rinse adhesive systems (groups 1 and 2) had significantly lower hardness and elastic modulus than those for region 1 (control enamel). The mechanical properties of enamel at least 6 μm from the external surface were influenced by demineralization, and this agrees with the results from previous studies of the effects of conventional etching and self-etching on the

Table 3. Mean Values and Standard Deviation for Cross-sectional Elastic Modulus (GPa) of Enamel at Distances From the Surface^a

Depth, μm	Region 1		Region 2		Region 3		Region 4	
	Mean	SD	Mean	SD	Mean	SD	Mean	SD
Group 1 : Transbond XT								
1	87.33 ^a	11.66	39.32 ^b	16.19	36.98 ^b	12.20	36.79 ^b	13.21
6	89.89 ^a	10.60	67.80 ^b	13.45	69.57 ^b	14.28	64.80 ^b	20.93
11	93.08 ^a	6.84	86.49 ^b	7.44	86.87 ^b	6.46	85.99 ^{a,b}	10.28
16	94.81	6.75	89.99	7.28	90.76	8.51	91.26	7.59
21	97.22 ^a	8.74	93.93 ^{a,b}	5.73	94.77 ^{a,b}	7.90	90.17 ^b	6.05
26	96.77	9.35	94.66	5.86	95.52	5.42	95.37	6.95
31	97.96 ^a	6.55	92.82 ^b	3.58	95.44 ^{a,b}	7.82	93.82 ^{a,b}	5.03
36	97.64 ^a	5.86	92.24 ^b	5.98	95.73 ^{a,b}	6.72	93.85 ^{a,b}	7.94
41	96.07	9.14	95.70	5.49	95.75	7.10	96.22	6.32
46	96.76	9.94	92.76	5.30	98.04	8.44	92.24	4.90
Group 2 : Transbond Plus Color Change								
1	95.44 ^a	5.58	59.89 ^b	27.56	53.03 ^b	18.99	54.47 ^b	21.88
6	92.95 ^a	7.38	70.44 ^b	20.70	72.47 ^b	18.33	72.10 ^b	20.16
11	93.34 ^a	5.69	87.83 ^{a,b}	9.48	85.94 ^b	6.91	88.80 ^{a,b}	8.87
16	96.02 ^a	6.19	91.19 ^b	6.30	86.54 ^c	8.38	87.69 ^{b,c}	9.04
21	96.57	5.86	91.85	7.36	94.13	6.21	91.74	8.97
26	95.98	4.55	93.61	4.51	97.26	8.01	93.79	7.37
31	97.11	4.87	94.66	6.73	97.31	6.17	95.79	7.74
36	98.93 ^a	4.93	93.07 ^b	5.40	97.93 ^{a,b}	5.77	95.44 ^{a,b}	10.36
41	99.29	8.92	97.02	6.79	98.28	5.67	94.29	13.56
46	100.34	7.94	96.36	6.70	98.22	5.51	96.59	6.86
Group 3 : Transbond XT + SEP								
1	91.07 ^a	5.75	88.95 ^a	8.08	74.48 ^b	9.93	76.91 ^b	11.90
6	93.23 ^a	6.33	95.14 ^a	4.61	78.18 ^b	8.52	80.07 ^b	10.22
11	94.62 ^a	4.30	95.37 ^a	6.42	83.46 ^b	9.75	81.89 ^b	13.05
16	96.41 ^{a,b}	4.05	99.13 ^a	5.52	88.61 ^b	10.27	88.70 ^{b,c}	9.91
21	98.45 ^a	7.09	97.02 ^a	4.21	91.73 ^b	5.73	90.35 ^b	6.42
26	96.52	5.65	98.00	7.30	94.25	6.72	94.76	6.29
31	97.39 ^{a,b}	5.77	99.50 ^a	8.61	94.13 ^{a,b}	7.37	92.63 ^b	5.56
36	98.09	7.21	98.90	6.97	93.66	8.15	93.34	5.85
41	99.28 ^a	7.61	101.18 ^a	10.06	92.57 ^b	5.98	94.15 ^{a,b}	7.35
46	96.81	7.57	98.58	6.49	96.78	6.21	98.21	8.83
Group 4 : Transbond Plus Color Change + SEP								
1	89.63	7.63	89.41	9.56	83.87	8.19	84.17	7.68
6	94.02	11.17	95.51	11.01	92.47	8.90	88.86	7.94
11	94.60	9.33	91.04	8.60	92.71	10.19	91.75	5.91
16	96.06	8.05	95.88	11.28	95.22	10.03	94.25	10.27
21	95.58	8.92	92.69	7.15	93.43	6.59	91.74	7.65
26	97.56	15.16	94.24	8.19	93.56	7.90	95.68	7.83
31	97.43	10.28	96.14	8.24	96.03	5.76	93.49	5.30
36	97.22	9.73	93.00	7.10	94.94	9.82	95.24	6.31
41	101.03	12.25	94.43	5.00	93.04	6.50	95.86	7.30
46	99.62	11.64	95.13	8.58	96.40	7.82	98.66	8.48
Group 5 : Beauty Ortho Bond								
1	91.97	10.35	87.10	8.07	85.00	11.87	82.91	11.26
6	91.41	10.24	93.52	8.43	90.77	10.09	91.04	9.04
11	93.45	11.58	95.59	9.32	95.10	9.26	92.29	9.76
16	93.73	9.76	94.15	8.88	95.27	10.30	96.59	11.37
21	93.45	10.04	95.55	8.10	93.44	8.41	99.37	7.94
26	93.91	8.66	97.32	8.51	97.34	10.35	96.58	7.60
31	97.05	9.39	93.25	7.39	97.83	8.10	100.43	10.63
36	97.75	6.86	98.03	9.21	97.48	7.94	97.66	6.74
41	97.80	6.27	96.02	6.65	97.16	7.64	98.76	6.21
46	97.37	6.04	96.59	7.05	97.02	5.94	99.72	10.22
Group 6 : Fuji Ortho LC Automix								
1	91.49 ^a	11.35	88.88 ^{a,b}	12.20	80.98 ^b	10.28	82.31 ^{a,b}	9.58
6	96.26	12.27	92.32	8.46	90.70	7.83	93.38	12.11

Table 3. Continued

Depth, μm	Region 1		Region 2		Region 3		Region 4	
	Mean	SD	Mean	SD	Mean	SD	Mean	SD
11	94.76	10.17	92.34	8.30	94.77	12.47	92.94	7.99
16	94.59	8.93	94.39	8.58	93.02	10.59	94.47	7.63
21	94.56	6.59	91.79	9.58	94.52	10.46	92.36	7.28
26	95.54	6.61	93.58	7.30	95.37	7.98	93.61	7.60
31	95.28	7.52	94.80	7.58	96.41	9.28	93.90	10.56
36	98.90	8.38	94.07	3.95	92.63	7.09	93.79	6.20
41	97.20	10.43	94.59	7.79	94.62	6.94	95.71	10.08
46	101.22	15.37	97.44	5.59	95.21	9.83	95.20	11.22

^a Region 1, control; region 2, under bonding adhesive; region 3, 50 μm from the bracket edge; region 4, 100 μm from the bracket edge. SEP indicates self-etching primer. Identical letters indicate that mean values were not significantly different ($P < .0125$).

mechanical properties (hardness and elastic modulus) of the enamel surface region.^{3,8} The location at 1 μm from the enamel surface for regions 2, 3, and 4 of groups 1 and 2 had similar values of hardness and elastic modulus, which was also found at the 6 μm location, and additional deterioration of the enamel structure after etching with 35% phosphoric acid was not observed. These results might suggest that Transbond Plus Color Change with fluoride-releasing ability was insufficient for remineralization of etched enamel by 35% phosphoric acid around brackets and that additional application of fluoride-containing products such as toothpaste, bonding materials, and solutions should be needed for the remineralization.

Over the past decade, progress has been made in bracket bonding to enamel with RMGICs and SEPs,² and their lower etching ability might minimize the potential for iatrogenic damage to enamel. Region 2 for groups 3, 4, 5, and 6 had higher values for hardness and elastic modulus than for groups 1 and 2, suggesting that the surface treatment of SEP and RMGIC adhesive systems causes minimal enamel demineralization; this agrees with the observation from a previous study.³ A previous study¹⁶ found that the etching patterns of aprismatic enamel were dependent on the aggressiveness of the acid, but there was no correlation between the aggressiveness and bond strength to intact enamel. Another study¹⁷ measured pH values of phosphoric acid and orthodontic SEPs, and the values were 1.39 for 35% phosphoric acid, 1.85 for Transbond SEP, and 2.20 for Beauty Ortho Bond SEP; both SEPs with relatively less acidic pH values had a mild etching effect for intact enamel. In this study, greater demineralization by 35% phosphoric acid compared to both SEPs was due to lower pH value and longer application time.

For group 3, the hardness up to 16 μm in depth and elastic modulus up to 21 μm in depth from enamel surface for regions 3 and 4 generally had significantly lower values than those for regions 1 and 2 (although there were no significant differences for region 1 at 16 μm depth in group 3), indicating that additional deterioration of the enamel structure by immersion in demineralization solution occurred after priming with SEPs. On the other hand, noteworthy deterioration of the enamel structure by immersion in the demineralization solution did not occur for groups 4 and 5. These results may confirm that the fluoride-releasing ability for groups 4 and 5 prevent enamel demineralization around orthodontic brackets.

Glass ionomer cements have high levels of fluoride-releasing ability, although they show poor bracket bond strength and greater bond failure rates than composite resin adhesives.^{6,18} To obtain adequate bond strength and provide greater fluoride release, RMGICs have been developed. The recently introduced RMGIC, Fuji Ortho LC Automix (GC America), in easy-to-use Automix paste packages, was selected for the present study. While Fuji Ortho LC Automix is considered to be an improved product based on Fuji Ortho LC adhesive (GC America), neither the fluoride release behavior nor the bracket bond strength have been reported for this new adhesive. Previous studies have reported that RMGICs showed more resistance to demineralization of enamel than other bracket bonding materials.^{1,6,7} However, SEP composite resin adhesive systems (groups 4 and 5) with fluoride-releasing ability showed an equivalent effect for preventive enamel demineralization in the present study.

Further materials-oriented studies are needed to provide insight about the relationships between bonding materials and enamel structures during orthodontic treatment. Such investigations may also provide new ideas for novel design of adhesive materials.

CONCLUSIONS

Under the conditions of this study, the following conclusions can be drawn:

- The hardness and elastic modulus of the enamel surface region is decreased by enamel etching with 35% phosphoric acid or priming with SEPs. The influence is minimal for SEP adhesive, and SEP adhesives or RMGIC adhesive systems with fluoride-releasing ability may increase the hardness and elastic modulus of the enamel surface region by remineralization of the enamel.
- The fluoride-releasing ability of the composite resin adhesive and RMGIC adhesive prevents enamel

demineralization and the deterioration in mechanical properties of enamel around brackets.

- Nanoindentation is a useful method for investigating mechanical properties in small regions of enamel, and these properties should be relevant to demineralization.

REFERENCES

1. Sudjalim TR, Woods MG, Manton DJ, Reynolds EC. Prevention of demineralization around orthodontic brackets in vitro. *Am J Orthod Dentofacial Orthop.* 2007;131:705.e1–e9.
2. Fjeld M, Øgaard B. Scanning electron microscopic evaluation of enamel surfaces exposed to 3 orthodontic bonding systems. *Am J Orthod Dentofacial Orthop.* 2006;130:575–581.
3. Iijima M, Muguruma T, Brantley WA, Ito S, Yuasa T, Saito T, Mizoguchi I. Effect of bracket bonding on nanomechanical properties of enamel. *Am J Orthod Dentofacial Orthop.* 2010;138:735–740.
4. Hu W, Zhou YH, Wang Q, Fu MK, Volpe A, Devizio W, Petrone M, Zhang YP. Effects of fluoride toothpaste on etched enamel of orthodontic patients. *Chin J Dent Res.* 1999;2:79–83.
5. Hu W, Featherstone JD. Prevention of enamel demineralization: an in-vitro study using light-cured filled sealant. *Am J Orthod Dentofacial Orthop.* 2005;128:592–600.
6. Paschos E, Kleinschrodt T, Clementino-Luedemann T, Huth KC, Hickel R, Kunzelmann KH, Rudzki-Janson I. Effect of different bonding agents on prevention of enamel demineralization around orthodontic brackets. *Am J Orthod Dentofacial Orthop.* 2009;135:603–612.
7. Gorton J, Featherstone JD. In vivo inhibition of demineralization around orthodontic brackets. *Am J Orthod Dentofacial Orthop.* 2003;123:10–14.
8. Pascotto RC, Navarro MF, Capelozza Filho L, Cury JA. In vivo effect of a resin-modified glass ionomer cement on enamel demineralization around orthodontic brackets. *Am J Orthod Dentofacial Orthop.* 2004;125:36–41.
9. Ghiz MA, Ngan P, Kao E, Martin C, Gunel E. Effects of sealant and self-etching primer on enamel decalcification. Part II: an in-vivo study. *Am J Orthod Dentofacial Orthop.* 2009;135:206–213.
10. Featherstone JD, ten Cate JM, Shariati M, Arends J. Comparison of artificial caries-like lesions by quantitative microradiography and microhardness profiles. *Caries Res.* 1983;17:385–391.
11. Oliver WC, Pharr GM. An improved technique for determining hardness and elastic modulus using load and displacement sensing indentation experiments. *J Mater Res.* 1992;7:1564–1583.
12. He LH, Swain MV. Influence of environment on the mechanical behaviour of mature human enamel. *Biomaterials.* 2007;28:4512–4520.
13. Alcock JP, Barbour ME, Sandy JR, Ireland AJ. Nanoindentation of orthodontic archwires: the effect of decontamination and clinical use on hardness, elastic modulus and surface roughness. *Dent Mater.* 2009;25:1039–1043.
14. Iijima M, Muguruma T, Brantley WA, Yuasa T, Uechi J, Mizoguchi I. The effect of filler on the grindability of composite resin adhesive. *Am J Orthod Dentofacial Orthop.* 2010;138:420–426.
15. ISO 14577-1. *Metallic materials—instrumented indentation test for hardness and materials parameters—Part 1: test method. 1st ed.* International Organization for Standardization; 2002.

16. Pashley DH, Tay FR. Aggressiveness of contemporary self-etching adhesives. Part II: etching effects on unground enamel. *Dent Mater*. 2001;17:430–444.
17. Iijima M, Ito S, Yuasa T, Muguruma T, Saito T, Mizoguchi I. Bond strength comparison and scanning microscopic evaluation of three orthodontic bonding systems. *Dent Mater J*. 2008;27:392–399.
18. Basdra EK, Huber H, Komposch G. Fluoride released from orthodontic bonding agents the enamel demineralization in vitro. *Am J Orthod Dentofacial Orthop*. 1996;109:466–472.